SIGMOIDOSCOPE
WITH LIGHT GUIDE, LENS AND INSUFFLATOR

Product Description

Description	Model Size	Product Code
SIGMOLUX*	One Size	701000

SIGMOLUX™ is a single use, fully disposable and self-illuminating rigid Sigmoidoscope. These devices are used during either rectal examination or rectal surgery, for use in examination, diagnosis and treatment of the large intestine. SIGMOLUX™ is supplied with an attached insufflator, a pre-loaded introducer and a hinged lens. The insufflator incorporates a hand squeeze pump with a one-way valve and reservoir (bladder) and is intended to inflate the colon during procedures according to standard clinical practice. All SIGMOLUX™ devices are shipped ready for use. Should the light output become obscured due to debris, the end of the light guide can be cleared by re-inserting the introducer. Re-use of this device represents an unnecessary and avoidable risk of cross infection to the patient.

Application

The G+N Medical SIGMOLUX™ and associated Illuminator are intended to provide a means of endoscopic access to the patient's large intestine and to supply effective illumination, thus permitting direct visual examination by a suitably qualified and experienced Healthcare Professional.

The devices key functions are:

1. To enable trained Healthcare Professionals, to effectively view patients' large intestine to perform invasive investigational procedures.
2. Reduce the complexity and increase the ease of use by having a light source directly attached to the device, enabling the large intestine to be viewed by the operator. The internal light source illuminates the area without the need for extra equipment.

Intended Use

The G+N Medical SIGMOLUX™ is intended to provide a means of rigid endoscopic access to the patient's rectum and large intestine. SIGMOLUX™ supplies effective illumination thus permitting direct visual examination by qualified medical personnel within a professional health facility. An insufflator is attached to assist in this procedure accordingly.

Indication and Patient Group

SIGMOLUX™ is indicated for use when a suitably qualified and experienced Healthcare Professional is required to perform a direct examination of the large intestine of a non-contraindicated adult

Not For Paediatric Use

Instructions for Use

1. Examine the pack before use. Remove the device from its packaging by separating the peel pouch and removing the product. Discard the package. The pack shall not be issued nor the Sigmoidoscope used if the wrapper or seal is broken. DO NOT USE IF PACKAGING OR SIGMOIDOSCOPE IS DAMAGED.
2. Check that all parts are correctly attached before usage
3. Pull out the blue tab in the light pack located in the device handle to activate the light source, before insertion into the patient. Maximum illumination intensity is maintained for a minimum of 30 minutes. The light pack will operate for a considerably longer period, with gradually diminishing intensity.
4. Lubricate SIGMOLUX™ with a suitable lubricant starting from the tip of the introducer before use on the patient.
5. It is advised and may be necessary to empty the rectum with an enema prior to performing the procedure. The patient should be informed when the Sigmoidoscope is to be introduced and removed.
6. Hold SIGMOLUX™ by the handle, with your thumb press against the introducer to avoid introducer shifting back and gently insert the instrument.
7. Once SIGMOLUX™ is placed in position, remove the introducer, close the lens and lock using the lens latch.
8. The instrument will be airtight when the lens is successfully closed.
9. Start insufflating and move SIGMOLUX™ towards the patient with lateral movements.
10. The lens window can be opened and closed at any time
11. On completion of the investigational procedure, press the button on the insufflator to release air from the colon. Remove SIGMOLUX™ from the patient with a gentle, lateral movement.
12. Dispose of the device in an appropriate way as described in the disposal section. **IMPORTANT** The lens can be used as protection against any unwanted material falling inside the Sigmoidoscope. If the examination is not being conducted and the Sigmoidoscope is still inserted, it is highly recommended to close the window. **IMPORTANT** Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Contraindications

Acute peritonitis	Acute Severe Diverticulitis	Diverticular Abscess
Colonic Necrosis	Recent Colonic Surgery	Anal Fissure
Fulminant Colitis	Toxic Megacolon	

	Do not re-use		Do not use if package is damaged
	Catalogue number		Keep dry
	Batch number		Shock protection (Type BF Applied Part)
	Use by date		Federal or EU Law restricts this device to be sold by or on the order of a licensed Health Practitioner
	Date of manufacture		Warning/caution consult accompanying documents
	Direct electrical current at voltage shown		Consult instructions for use/consult operating instructions
	Not made with natural rubber latex		Non-sterile
	Do not look directly into the light source		Manufacturer and address of the manufacturer
	Fragile – handle package carefully		Number of items per pack
	No DEHP		IP rating
	Medical Device		Medical Device in conformity with Regulation (EU) 2017/745

GMDN 15058	Global Medical Device Nomenclature
------------	------------------------------------

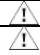
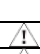
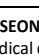
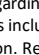
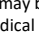
Disposal

Discard SIGMOLUX™ with the light pack still attached, as "clinical waste" according to hospital, or local protocols. Note: If required the light pack can be removed by firmly pulling the entire pack downwards. Do be careful as it may well release suddenly. Dispose of the light pack responsibly.

Warnings

Before each use carefully examine the instrument to ensure that there are no visible scratches, chips, voids, etc. These could cause the instrument to lose its stability potentially compromising safety. The patient should be informed when the Sigmoidoscope is to be introduced and removed.

Adjustments may be necessary to provide improved comfort.

WARNING		Inspect each device and packaging for damage prior to use
WARNING		For hygienic reasons always wear surgical gloves when handling the device before and after use. The device will come into contact with biological fluids which are potentially infectious.
WARNING		When in transit or storage, device may be subject to damage beyond the control of the manufacturer or supplier.
WARNING		Never use device with laser equipment.
WARNING		Treat used devices as biohazardous infectious material. Dispose of used devices in suitable disposal unit or in accordance with local regulations.
WARNING		Device is not compatible with any sterilization method.
WARNING		No modification of this equipment is allowed.
WARNING		Never use an ultrasonic cleaner for the device, the device should not be cleaned or sterilised in any way as it provided in clean hygienic form for single use only.

FOR SINGLEUSEONLY:

The G&N Medical device is designed for single use only; do not reuse device. G&N does not have data regarding reuse of this device. Reuse may cause device failure or procedural complications including device damage, compromised device biocompatibility, and device contamination. Reuse may result in infection, serious injury, or patient death. After use, this product may be a potential biohazard. Handle and dispose of in accordance with accepted medical practice and applicable laws and regulations.

Cautions:

For use by trained healthcare personnel only

Federal (USA), UK and European Union Law restricts this device to sale, distribution and use by or on the order of a physician.

Storage

Store devices in their original packaging at room temperature and normal air humidity.

CS Lifesciences Europe Limited 3 Inns Quay, Dublin 7, Ireland

110001/3
05/2021

Griffiths and Nielsen Ltd.

Maydwell Avenue, Stane Street, Slinfold, Horsham, RH13 0GN, UK

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