


SIZING CHART

(23-32mmHg)



1

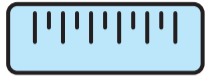

Measure ankle circumference in centimetres (cm)
Two-fingers width above ankle bone

2



Check ankle measurement chart for size

3




Observe limb proportions
Is the calf in proportion to the ankle?

4




Are you fitting below knee or thigh length?

YES

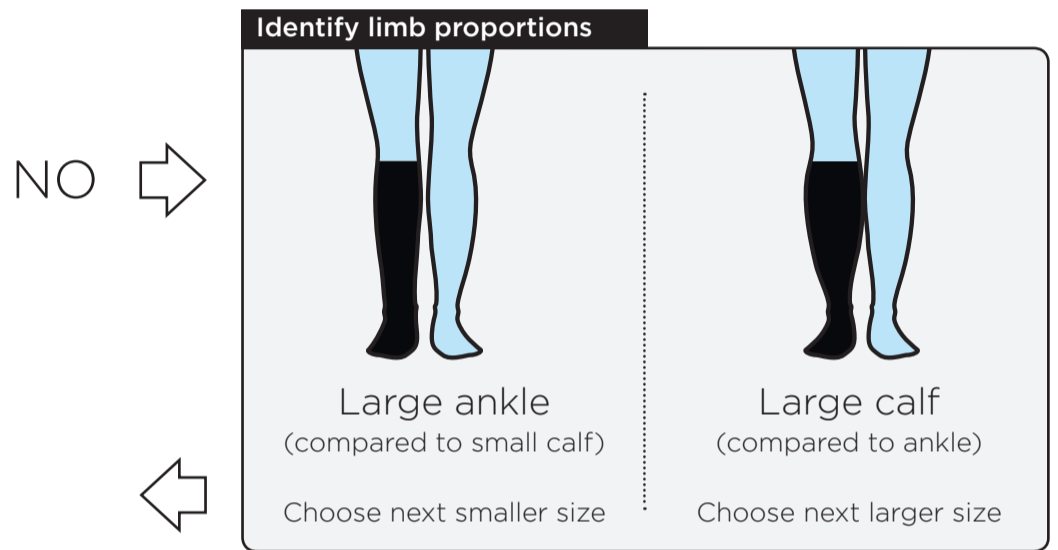


Refer to packaging for application instructions and contra-indications

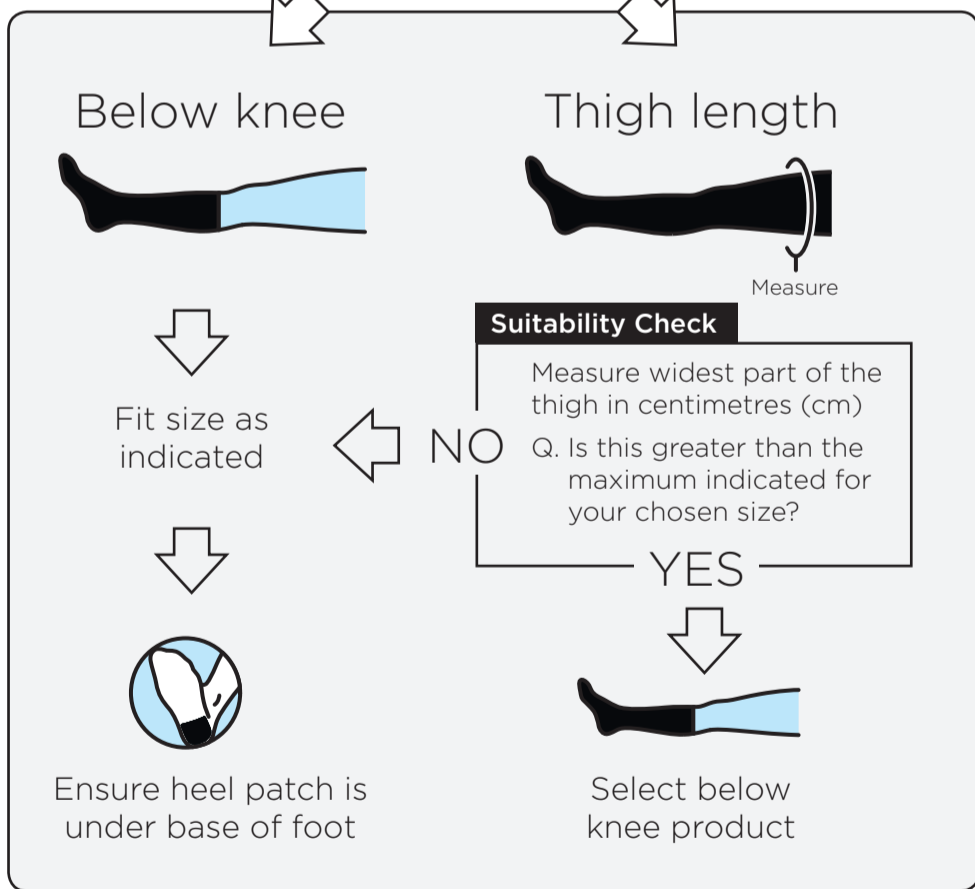

FITLEGS® 2 (23-32mmHg)

Suitability Check

	Ankle measurement	Size	BELOW KNEE Order Code	THIGH Order Code	Measurement based on height & weight	TIGHTS Order Code	Maximum Thigh Circumference
Black	19 - 23cm	Small	457030	457034	See table below for full sizing for Fitlegs Class 2 tights range	457046	65cm
	23 - 26cm	Medium	457031	457035		457047	70cm
	26 - 29cm	Large	457032	457036		457048	78cm
	29 - 33cm	Ex-Large	457033	457037		457049	82cm
Beige	19 - 23cm	Small	457016	457021	See table below for full sizing for Fitlegs Class 2 tights range	457050	65cm
	23 - 26cm	Medium	457017	457022		457051	70cm
	26 - 29cm	Large	457018	457023		457052	78cm
	29 - 33cm	Ex-Large	457019	457024		457053	82cm
White	19 - 23cm	Small	457038	457042	See table below for full sizing for Fitlegs Class 2 tights range		65cm
	23 - 26cm	Medium	457039	457043			70cm
	26 - 29cm	Large	457039	457044			78cm
	29 - 33cm	Ex-Large	457040	457045			82cm

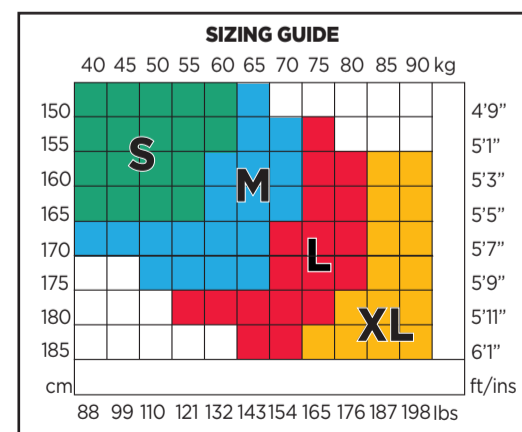

5



Record chosen size and document in notes
Re-check measurements and fitting regularly